




Bipolar Eyelash Pencil



Unique Micro-fine 23 gauge Bipolar Pencil for eyelash removal and minor in-office procedures.

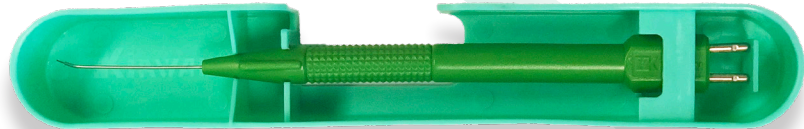

14-5137 Disposable 23 gauge bipolar pencil, bent, with 3mm extended tip. This unique bipolar pencil is 23 gauge with a 3mm 30 gauge extended tip. The knurled handle offers comfort and precise control while the 45 degree bent shaft increases access to target tissues even while the patient is in a slit lamp. The 14-5137 bipolar pencil is safer to use, is associated with less post operative pain and promotes quicker recovery times. It is very low cost, easy to set up and use as well as compatible with any standard bipolar generator and standard bipolar cord.

Other Kirwan products to use with the 14-5137 disposable 23 gauge bipolar pencil:

- 26-2500 AURA 110/220V 20 watt bipolar generator
- 10-5000 Reusable 12' TPR bipolar cord, autoclavable
- Package purchases with generator, pencils and cords are available.

